

emergency

# Silicone Resuscitators

*Clear Ventilation*



# Silicone Resuscitators

## *Clear Ventilation – Increased Safety for Your Patients*

### **Flexible**

The rotating patient connection lends the bag-valve-mask greater flexibility and its structured surface gives the user a good grip. Because it is available in three different sizes, you'll always have the right silicone resuscitator ready to treat adults, children and babies. You also have the choice from a range of accessories, such as a demand valve or an oxygen reservoir for additional oxygenation.

### **Safe**

Thanks to the transparent mask, you can always keep an eye on your patient. If the patient vomits, for example, you can see the problem immediately and react without delay.

### **Clean**

In the transparent bellows and mask it's easier for you to see traces of contamination. That makes hygienic conditioning before the next use easier too.



Safety valve



### **Article No.**

- |   |                 |
|---|-----------------|
| <b>1 Silicone bag-valve-mask for babies/toddlers</b>  | <b>WM 11101</b> |
| 240 ml, including oxygen reservoir 600 ml and three silicone ventilation masks, sizes 0, 1, 2 |                 |
| <b>2 Silicone bag-valve-mask for children</b>   | <b>WM 11102</b> |
| 500 ml, including oxygen reservoir 2600 ml and silicone ventilation mask, size 3              |                 |
| <b>3 Silicone bag-valve-mask for adults</b>   | <b>WM 11103</b> |
| 1600 ml, including oxygen reservoir 2600 ml and silicone ventilation mask, size 5             |                 |



Adapter fits oxygen reservoir  
or OXYMAND demand valve

#### Advantages and benefits:

- Safety assured by integrated pressure limit valve
- Simple hygienic conditioning: completely autoclavable
- Transparent: for effective patient monitoring and fast detection of contamination
- Several connection options, such as to oxygen reservoir or demand valve
- Latex-free: all components are of silicone or other latex-free materials
- Always fits: the silicone bag-valve-mask is available in three different sizes
- Rugged construction: suitable for the toughest operating conditions

1



2



3



## Silicone resuscitators



### Accessories and replacement parts for silicone resuscitators

■ Silicone bag-valve-mask (Rendell-Baker)			
1 Size 0 for babies	WM 5060	6 PEEP valve, infinitely adjustable from 0-10 mbar	WM 3215
2 Size 1 for toddlers up to one year	WM 5061	■ Respiration hose, steramid, Ø 6/9, length 2000 mm	WM 1411
3 Size 2 for children one to three years old	WM 5062	■ Patient valve	WM 11104
4 Size 3 for children three to 12 years old	WM 5063	■ Reservoir valve	WM 11105
5 OXYMAND demand valve with 1.5-meter pressure hose and adapter for Weinmann's silicone bag-valve-mask		■ Oxygen reservoir bag 600 ml	WM 11106
■ With G 3/8" straight connection and swivel nut for Weinmann pressure reducer	WM 22190	■ Oxygen reservoir bag 2600 ml	WM 11107
■ With Walther connector	WM 22167	■ Intake valve, complete for silicone bag-valve-mask, all three sizes	WM 11120
■ With DIN 13260 plug for central gas facility	WM 22168		
■ Adapter for OXYMAND demand valve	WM 22187		


### Technical data silicone resuscitators



<b>Product class as per directive 93/42/EEC:</b>	IIa	<b>Dimensions:</b>	
<b>Area of application:</b>		■ Babies (WM 11101):	8 x 13 x 48 cm
■ Babies/toddlers (WM 11101):	weight of 0 - 10 kg	■ Children (WM 11102):	9 x 15 x 62 cm
■ Children (WM 11102):	weight of 10 - 40 kg	■ Adults (WM 11103):	12 x 18 x 69 cm
■ Adults (WM 11103):	weight > 40 kg	<b>Weight:</b>	
<b>Pressure limits:</b>	40±5 cm H <sub>2</sub> O or unlimited	■ Babies (WM 11101):	300 g
<b>Gas inlet connector:</b>	Ø 6 mm outer diameter	■ Children (WM 11102):	370 g
<b>Patient connection:</b>	22/15 mm (ANSI / ISO)	■ Adults (WM 11103):	530 g
<b>Materials:</b>	Silicone rubber, polycarbonate, Polyvinyl chloride, stainless steel, polysulfone All materials latex-free	<b>Volumes:</b>	
		■ Babies (WM 11101):	240 ml
		■ Children (WM 11102):	500 ml
		■ Adults (WM 11103):	1600 ml

### Simply Professional

You can rely on WEINMANN Emergency. We set the standards for mobile solutions in patient ventilation. Our 130 years of experience and close association with professional users in emergency medical services, armies/public authorities and hospitals have led to the development of intelligent and perfectly coordinated technologies. In conjunction with the professionals in emergency, transport and disaster medicine, we pursue our shared goal of saving human lives.

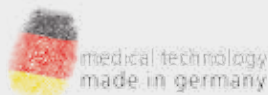

weinmann-emergency.de | weinmann-emergency.fr | weinmann-emergency.ru | weinmann-emergency.cn

Would you like to know more?  
Call us on  
**+49 40 88 18 96-120**

